


ANIMAL HYPOXIC GLOVE BOX



01 ACCURACY

- Precise O₂ control
 - ✓ Range of O₂ Control: 0.5% - 20.9%, in 0.1% increments
- Precise temperature control (Optional)
 - ✓ Range of Temperature Control: 5°C above ambient up to 45°C, in 0.1°C increments
- Advanced Airway Design
 - ✓ Ensures a high level of homogeneity in the workstation's gas distribution

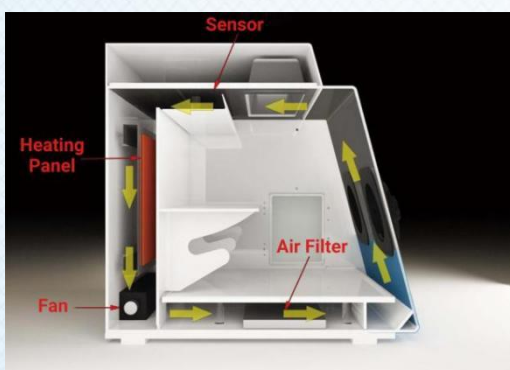

02 CONVENIENCE

- Easy Access Operation System
 - ✓ The sleeves are made of a soft, skin-friendly material and can access the workstation easily and quickly without vacuum

- **Optional Direct-access operation system**

- ✓ The sleeves are made of a soft, skin-friendly material and an efficient vacuum pump direct access to the workstation only 8-12 seconds.
- ✓ Can be operated directly with lab gloves.
- ✓ Glove port caps can be "parked" in special holders and do not take up valuable working space.



Direct-access Operation System



- Interlock can be quickly purged by simply pressing the yellow button

- An optional built-in air pump and filter are used to provide the oxygen supply, avoiding the need for a compressed air cylinder

03 CUSTOMISED PRODUCT

